


Createch for Straumann®

Customized prosthetic solutions.
Individual excellence and efficiency
on demand.



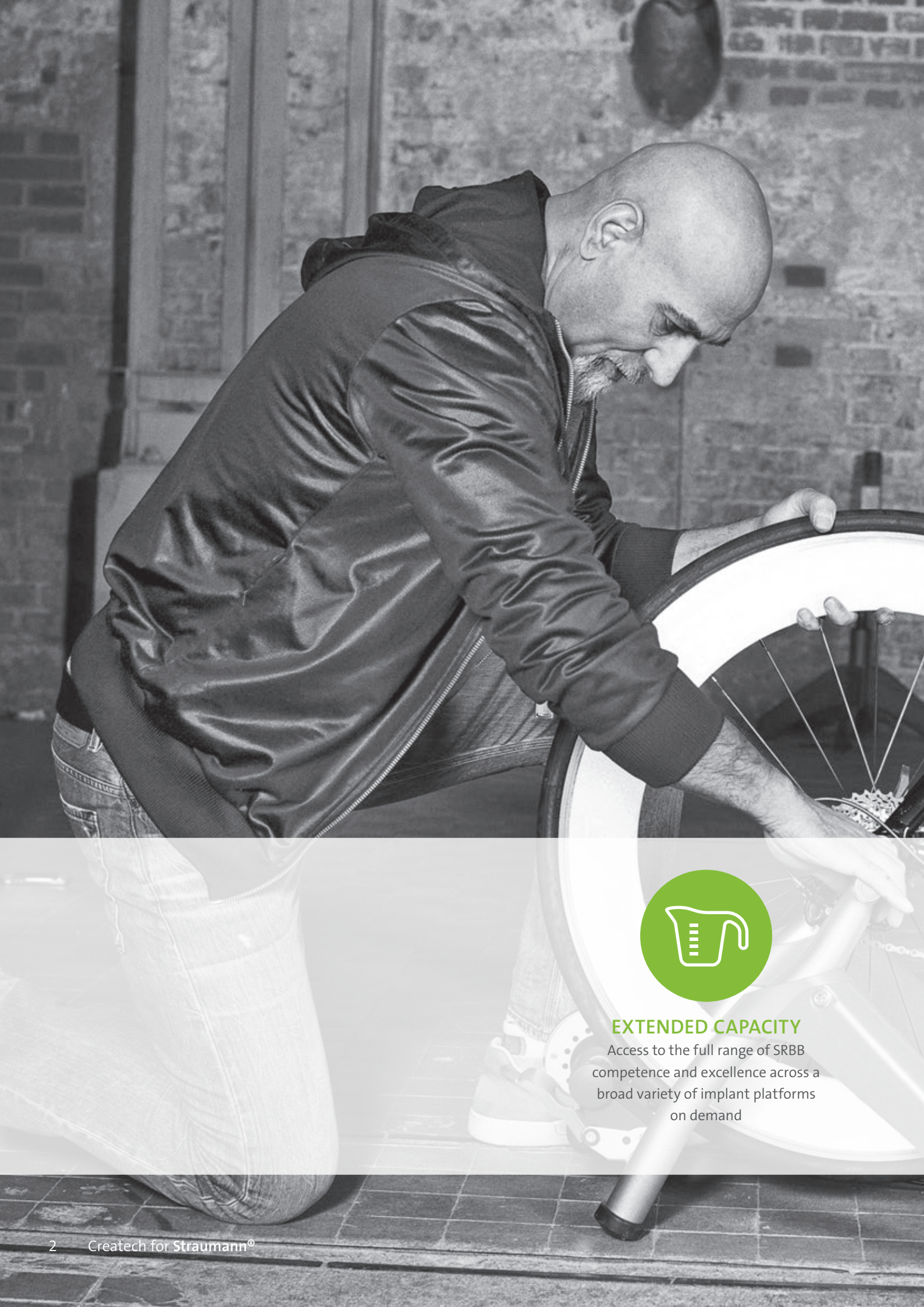

EXTENDED CAPACITY

Access to the full range of SRBB competence and excellence across a broad variety of implant platforms on demand

Customized prosthetic solutions. Individual excellence and efficiency on demand.

Customized screw-retained bars and bridges (SRBB) for complex restorations requiring a high degree of individualization are the equivalent of the “Champions League” for dental technicians and professional restorative experts – often difficult to be achieved with one’s own resources, with the necessary precision or in a given timeframe.

The Createch for Straumann® product range gives demanding customers extended capacity, competences and flexibility on demand. Createch is a Straumann Group certified partner offering customized prosthetic premium solutions for almost all implant systems available in Europe.



CONVENTIONAL AND DIGITAL

Unparalleled service and support on demand through Createch’s modular Excellence and Service Line that lets you choose your most convenient workflow



CUSTOMIZED EXCELLENCE

Tailor-made solutions for maximum flexibility: High degree of individualization to best meet your patients’ needs, including angulated screw channels



PREMIUM PRECISION

Consistently high quality and a highly precise fit for your restorations



Extended capacity: Createch Medical and Straumann

ACCESS TO THE FULL RANGE OF SRBB COMPETENCE AND EXCELLENCE ACROSS A BROAD VARIETY OF IMPLANT PLATFORMS ON DEMAND

Createch Medical is a Straumann Group member and stands for high-quality screw-retained bars and bridges for fixed and removable prosthetics. As a pioneer in its field, Createch Medical was founded in 2004 with the goal of meeting an existing need in the market for CAD/CAM dental frameworks.

Over the years the company has developed a comprehensive range of prosthetic solutions that facilitate the work of both dental technicians and professional restorative experts. At the same time the high degree of customization offered by Createch products enables dental professionals to provide tailor-made solutions for their patients.

The uniqueness of its portfolio offering, combined with its unrivalled quality and outstanding service and support, have made Createch a leader in its field.

With Createch for Straumann®, Createch supplements Straumann's portfolio offering of screw-retained bars and bridges with complex restorations requiring a high degree of individualization. Furthermore, dental professionals can now also benefit from Createch's compatible frameworks that can be used with over 30 implant brands, and over 200 implant connections.

HIGHEST QUALITY AND PRECISION

Createch fuses the precision of the aeronautical technology of an established industrial group with the creativity of implantology experts and laboratory technicians in order to provide the best possible outcome for its customers. Ultra-high precision measuring devices and state-of-the-art milling equipment, combined with strict quality standards in the manufacturing process, from scan to final restoration, ensure consistently high quality and accuracy of fit for your Createch frameworks.



Conventional and digital: Unparalleled service and support

Createch for Straumann® offers you the possibility to send your cases in your most convenient and efficient way, regardless of your digital capabilities. This will free up valuable time and resources in your lab, enabling you to optimize your productive time.

If you prefer the conventional workflow, choose Createch's Excellence Line and let Createch do the scan and design for you. Or if you prefer to do the scan and design with your own scanner, simply send us your design data.

Guarantees:

IN CHOOSING CREATECH FOR STRAUMANN® YOU ARE OPTING FOR A PREMIUM SOLUTION OF THE HIGHEST QUALITY FOR YOUR PATIENT

Benefit from the combined Createch and Straumann guarantee

Createch restorations on Straumann implants comes with the original Straumann connection. Your benefit is the 10-year guarantee* offered by Createch for the prosthetic framework. Plus, when restoring Straumann implants with a Createch solution, you also enjoy the lifetime (plus) guarantee** on the Straumann implants.

Enjoy the Createch guarantee when restoring other implant systems

Createch offers a 10-year guarantee on prosthetic frameworks compatible for other implant systems. In addition, Createch reimburses the cost of the replacement implant if the original implant manufacturer refuses its guarantee due to combination with a Createch product.*



&



*for more info see the Createch terms and conditions.

**for more info see the Straumann terms and conditions.

ONE LOGIN, TOTAL FLEXIBILITY

Visit www.Createchmedical.com to:

- Create user ID
- Download Createch library
- Create cases
- Choose preferred workflow line
- Monitor your case history
- Order DHL pick-up

Option “Service Line”

Provides a series of services to laboratories with a scanner:

- Training and consultation for the laboratory
- Scanner optimization
- Incorporation of our Createch solutions
- High-precision milling

Steps:

- Scan and design the framework in your preferred CAD environment
- Make sure the correct scan bodies are used
- Create and upload the case
- Send in physical model for precision scanning
- Create DHL button for pick-up order
- Print delivery note
- Allow 5 days turnaround time (3 days processing* plus shipment time depending on shipment option chosen)

Option “Excellence Line”

For those who don't work digitally and want the full service:

- Measuring the implants
- Designing the structures
- Simulation using finite element calculation
- Milling using the most advanced technology

Steps:

- Create the case
- Send in physical mock-up or framework, including the model to be processed
- Create DHL button for pick-up order
- Print delivery note
- Allow 7 days turnaround time (5 days processing* plus shipment time depending on shipment option chosen)

*Processing starts as soon as Createch has received all parts and information needed





Customized excellence: Maximum flexibility

TAILOR-MADE SOLUTIONS FOR THE ULTIMATE IN FLEXIBILITY

Createch for Straumann® provides you with all the flexibility you need in order to optimally restore your prosthetic case with a unique tailor-made solution, no matter how simple or complex it may be. Not only can you choose whether to restore your case at implant or abutment level, you can also opt for angulated screw channels, up to 25° angulation, in order to bring the screw exit hole into a favorable position. Moreover, Createch for Straumann® also offers you the opportunity to solve cases either with Straumann implants or other implant systems.



Required torque in Newton centimeter (Ncm) for your Createch restoration

Bridge, implant level, straight channel	35 Ncm
Bridge, abutment level, straight channel	15 Ncm
Bridge, implant or abutment level, angulated channel	15 Ncm



In case you need scanbodies on non-Straumann platforms, contact your local Straumann representative for help.

Digitalize your case:

Choose from a list of available scanbodies when restoring Straumann implant cases



Tissue Level – implant level



Straumann® CARES® Mono Scanbody*

NNC	 048.173
RN	 048.168
WN	 048.169



Bone Level – implant level



Straumann® CARES® Mono Scanbody*

NC	 025.2915
RC	 025.4915











Bone Level – abutment level (on Straumann® Screw-retained Abutments)

Straumann® CARES® Mono Scanbody*

\varnothing 3.5 mm	 025.0000
\varnothing 4.6 mm	 025.0001

Premium precision: Consistently high quality and a highly precise fit for your restorations

	Hybrids			Telescopic solutions		
Createch Screw-retained restorations*						
	Hybrid	Abutment hybrid	Basic hybrid	Removable telescopic	Removable anatomic telescope	Telescopic hybrid
	Restoration on implant level	×	×	×	×	×
	Restoration on abutment level*	×	×	×	×	×
Material						
Cobalt chromium			×			
Titanium Grade 5	×	×	×	×	×	×

	Bars			
Createch Removable restorations				
	Createch Bar	Dolder® Bar	LOCATOR® Bar	Preci-Vertex® Bar
				
	Shape			
Restoration on implant level	×	×	×	×
Restoration on abutment level*	×	×	×	×
Material				
Cobalt chromium	×	×	×	×
Titanium Grade 5	×	×	×	×




× = available

Additional information

* restorations available for abutment level for the following platforms:
Straumann® Multibase, Straumann® Screw-retained Abutment, synOcta® 1.5

Createch may not be available in all countries



Screw-retained bridges		
		
Anatomic telescopic hybrid	Full-arch bridge	Screw-retained bridge with gum support
×	×	×
×	×	×
	×	×
×	×	×



International Headquarters

Institut Straumann AG

Peter Merian-Weg 12

CH-4002 Basel, Switzerland

Phone +41 (0)61 965 11 11

Fax +41 (0)61 965 11 01

www.Straumann.com

Dolder® is a registered trademark from Prof. Eugen Dolder, exclusively licenced to Cendres + Métaux Holding SA, Switzerland.
Preci-Vertex® is a registered trademark of Alphadent N.V., Belgium.

© Institut Straumann AG, 2017. All rights reserved.

Straumann® and/or other trademarks and logos from Straumann® mentioned herein are the trademarks or registered trademarks of Straumann Holding AG and/or its affiliates.