


STRAUMANN® CARES® DIGITAL SOLUTIONS

Straumann® CARES® Digital Solutions > Intraoral Scanner

Straumann® Virtuo Vivo™
Capture each note.

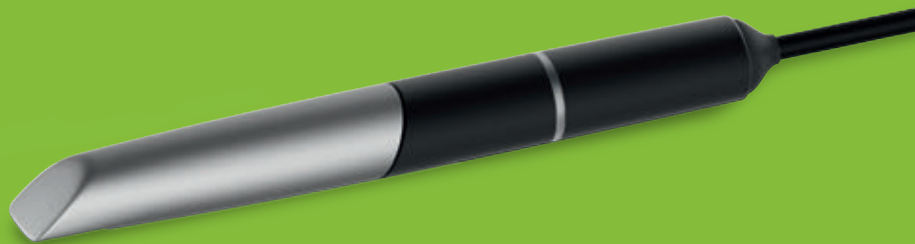

Virtuo Vivo™, the center of digital ecosystem

A COMPLETE INTEGRATED AND VALIDATED DIGITAL WORKFLOW: FROM SCAN TO MANUFACTURE

Accurate, powder free scanning in real life color is just a beginning of your digital journey. With Virtuo Vivo™ you have access to full digital ecosystem.

coDiagnostiX®

- More than implant planning software: a practice differentiator.
- coDiagnostiX™ offers a new groundbreaking guided surgery workflow which integrates both drill guide, prosthesis design and fabrication. Directly send the intraoral impressions from Virtuo Vivo™ to coDiagnostiX™ allows the planning and design of surgical guides in a perfect seamless way.



CARES® Visual Chairside One day dentistry

- Transfer the files for chairside design of restorations with CARES® Visual Chairside or with other 3rd party software.
- Color and accurate scans allow precise and user friendly design of restorations.
- Streamlined connection with C series chairside milling machine.

Connectivity

- With the open file format, you can send your digital impression directly to your lab partner via Straumann® CARES® Connect or direct STL or PLY export.
- Virtuo Vivo™ is an open system with validated workflows for ease of use.

Orthodontics

- Virtuo Vivo™ provides seamless integration with ClearCorrect.
- Send your files directly to ClearCorrect for orthodontics treatment with clear aligners.
- Since clear aligners are removable, and nearly invisible, many patients are actually excited to be treated with them.
- Virtuo Vivo™ direct integration to ClearCorrect facilitates submissions for clear aligners.

Design & production services

- You can also send the files for prosthetic design directly to Scan & Shape.
- Our team of highly skilled, digital dental technicians will prepare and produce the restoration in one of our centralized production centers.



Replace traditional dental impressions with accurate digital data

With Virtuo Vivo, digital impressions are acquired quickly, easily and in real color.



EASE OF HANDLING

Virtuo Vivo™ is designed to offer light and small handpiece. It is ergonomically designed to provide a pen-grip hold and comfort to the user. Additionally, the Virtuo Vivo™ comes with removable sleeves that eliminate patient turnaround time. The sleeves can withstand sterilization in an autoclave for up to 250 cycles.

UNSURPASSED SCANNING ACCESS

Multiscan Imaging™ captures data from many angles simultaneously. This 3D scanning technology, packs the power of two miniaturized 3D scanners into small hand-held intraoral scanner. Virtuo Vivo™ provides real life color scanning which allows for a better communication between the dentist and the patient and easier margin line marking.





GESTURE & VOICE CONTROL

Wave goodbye to traditional touchscreen monitors. Motion control technology allows touch-free manipulation of the screen imagery while wearing gloves.

The scanner is also equipped with a voice and a gesture control module. Using voice or gesture control instead of the touchscreen prevents cross-contamination by allowing to perform all scans without moving your hands away from the patient or the handpiece. Additionally you can use the air-mouse feature that enables movement of scan on the screen based on the tilt of the handpiece.

REAL-TIME DIGITAL RESULTS

A luminescent ring on the handpiece and audible signals allow the user to have confirmation that the data is being acquired. Data quality can be validated in real-time as the software creates the 3D model. The data can then be sent immediately to an external service provider.

PATIENT-FOCUSED DESIGN

Virtuo Vivo™ replaces conventional tray impressions, which are unpleasant for the patient, allowing you to discuss results and treatment options in real time.

CARES® CONNECT

CARES® Connect provides any member of the dental clinic team the ability to digitally send intraoral or impression scans to the laboratory of choice, and tracking tools will provide information on the status of the case.

We are here to help you take every step on your digital journey.

Straumann® CARES® Digital solutions can set the stage for the future of your practice, with plug and play simplicity.

With Virtuo Vivo™ you can expand your offering, delight patients and dial up your efficiency

- Best of all, you'll never rock alone with Straumann®
- Dedicated equipment specialists
- Thorough training and ongoing support





Contact your local
Straumann® representative
to schedule a **live demo**.



Straumann® Virtuo Vivo™

Straumann® Virtuo Vivo™

Rock your practice

- Scanning in real color: Increased field and depth of view
- Small and light handpiece
- Removable & autoclavable sleeves eliminate patient turn around time



490486/en/A/00 03/19

International Headquarters

Institut Straumann AG
Peter Merian-Weg 12
CH-4002 Basel, Switzerland
Phone +41 (0)61 965 11 11
Fax +41 (0)61 965 11 01
www.straumann.com