

*SMALL DIAMETER,  
GREAT ACHIEVEMENTS.*

**NEODENT<sup>®</sup> HELIX GM NARROW**





# Neodent® Helix GM Narrow

Small diameter, great achievements.

Bring reliability to your practice through the next generation of immediate esthetic solutions for reduced interdental spaces and bone availability.

## Ø 2.9

Patient expectations for immediate and esthetic treatment for missing teeth are constantly growing, even in cases with challenging indications such as reduced interdental spaces and reduced bone availability.

Implant therapy in such situations can be perceived as challenging for even highly experienced clinicians, due to implant strength concerns and lack of adequate prosthetic options.<sup>(1)</sup>

Finally, surgical procedure approach, both guided or non-guided handling can get arduous due to complex treatment procedures.

The Ø 2.9mm **Helix GM Narrow** provides an immediate, small diameter solution seeks to provide simplicity for treatment protocol – regardless of whether guided or non-guided techniques are used – confidence without compromising on strength, and flexibility for immediate esthetic outcomes in limited interdental spaces.

“

“With the new **Helix GM Narrow** system, I’m confident that I’m using a reliable long-term solution that allows for treatment of patients with reduced interdental space, with expected esthetic results. The narrow implant design features adequate primary stability that enables the use of immediate loading. The versatility of the hybrid surgical kit, as well as the comprehensive prosthetic portfolio, allows me to customize surgical plans seeking to provide a predictable and successful results for my patients.”



DR. GENINHO THOMÉ

”




**CONFIDENCE WITH A STABLE LONG-TERM IMPLANT FOUNDATION**

A unique combination based on proven concepts.



**FLEXIBILITY FOR IMMEDIATE ESTHETIC OUTCOMES**

Versatile portfolio designed for limited spaces.

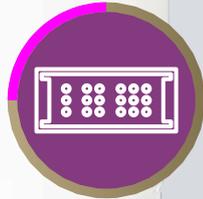

**SIMPLICITY FOR TREATMENT PROTOCOLS**

Intuitive and functional surgical cassette to best suit any chosen surgical procedure.





## Confidence with a **stable long-term implant foundation**


Implant therapy for demanding indications, such as reduced interdental spaces, can raise concerns regarding resistance and biomechanical behavior. Therefore, features of an implant-abutment interface are essential to provide successful long-term functional, stable, and esthetic results.

The Ø 2.9mm Helix features the strong and stable GM Narrow connection, designed with a unique combination based on proven concepts seeking to achieve long lasting results. A system produced out with the commercially pure titanium grade 4 offering treatment predictability through the Acqua hydrophilic surface.

### RELIABLE AND STRONG GM NARROW CONNECTION


#### 1 **16° Morse Taper connection**

The implant-abutment interface is a relevant aspect that could interfere on the success of patient's outcome. **Helix GM Narrow** is designed to deliver a tight fit for optimal connection sealing and offers strong mechanical resistance.

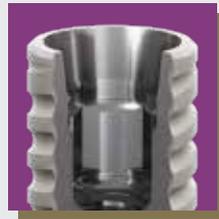

#### 2 **Internal hexagonal indexation**

The connection is designed with internal hexagonal indexation for precise abutment positioning, easy handling.

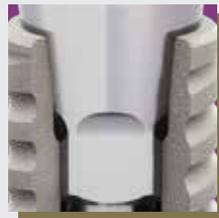


#### 3 **Platform switching**

The abutment design features a narrower diameter than the implant coronal area, which enables platform switching. <sup>(2-6)</sup>



#### 4 **Screw-retained interface**

The **Helix GM Narrow** features a morse taper screw-retained connection, which fits into the internal thread with precision seeking to provide a stable abutment connection.

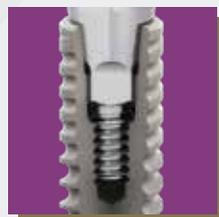

## COMMERCIALLY PURE AND MECHANICALLY STRONG TITANIUM GRADE 4

Beyond a versatile design allowing primary stability, the **Helix GM Narrow** is produced from the most commercially pure and mechanically strong titanium grade 4 (Ti Gr 4).

Static torsion tests have been conducted providing a greater performance and strongness of +12,7% than the former small diameter Neodent® system (Ti6Al4V-ELI).

**+ 12,7%**

### Static torsion test


New small diameter Neodent® system (Ti Gr 4)

Former small diameter Neodent® system (Ti6Al4V-ELI)


Font: Annex\_NoC Helix Narrow internal document.

## ACQUA HYDROPHILIC SURFACE'S AND TREATMENT PREDICTABILITY

The Neodent® Acqua hydrophilic surface is the next level of the highly successful S.L.A. surface. It was developed to reach expected results outcomes even in the most challenging patient cases, such as soft bone or immediate protocols. <sup>(7-10)</sup>

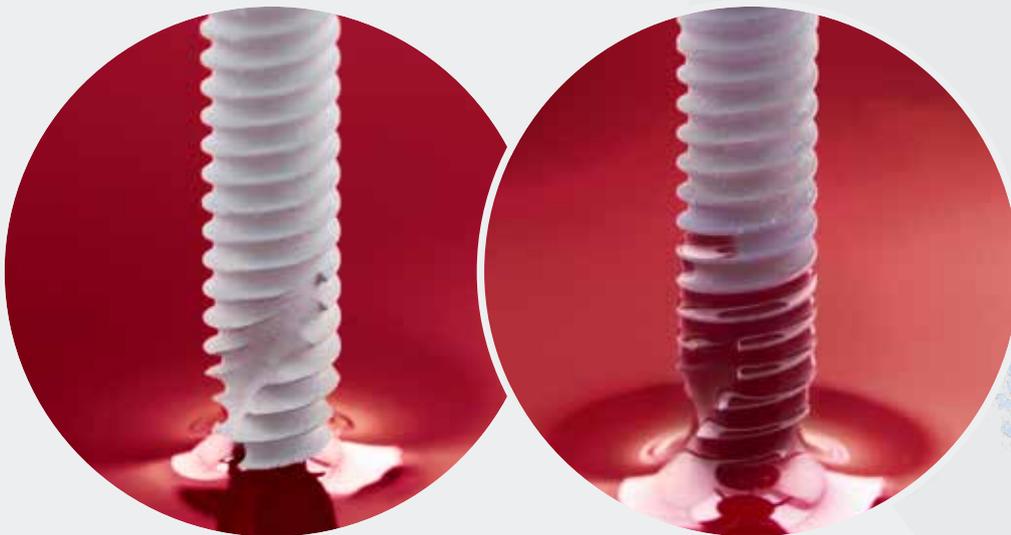

The Neodent® hydrophilic surface presents a smaller contact angle when in contact with liquids. This provides greater accessibility of organic fluids to the Acqua implant surface to help spur healing. <sup>(8)</sup>

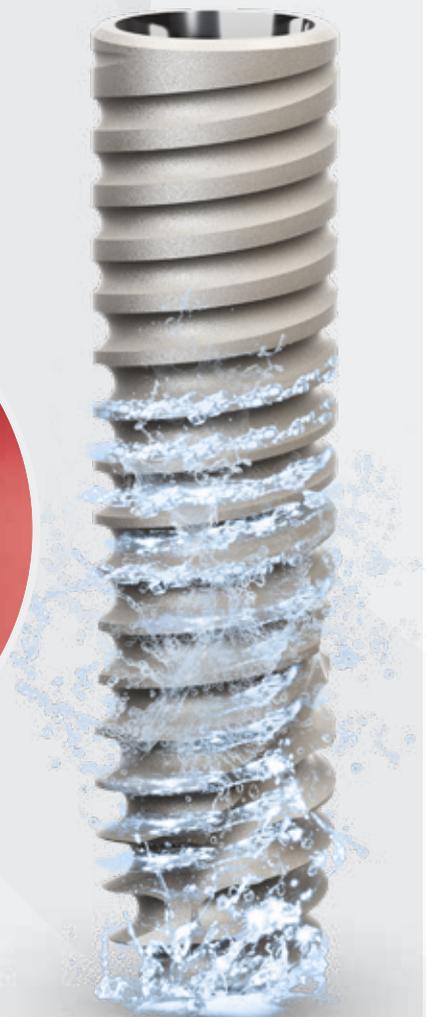



## Flexibility for immediate esthetic outcomes


Patients lacking bone availability in the esthetic zone or experiencing limited space between adjacent teeth, can make tooth replacement procedures challenging for implant clinicians. When coupled with a lack of adequate prosthetic options to correctly replace missing teeth, patient satisfaction declines, and practices can suffer.

The versatile **Neodent® Helix GM Narrow** system combines a  $\text{\O}2.9\text{mm}$  Helix implant, with a comprehensive prosthetic portfolio to restore cases in limited bone availability and interdental spaces, for immediate esthetic results.

### THE UNBEATABLE VERSATILITY OF HELIX


#### Dynamic progressive thread design

- Coronal: Double start threads with rounded root > compressing;
  - Apex: V-Shape > Self-cutting
- High primary stability.

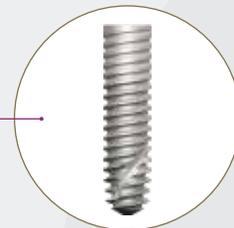


#### Tapered body design

- Coronal: Progressive tapered design;
  - Apex:  $12^\circ$
- Under-osteotomy for bone types 3 and 4.

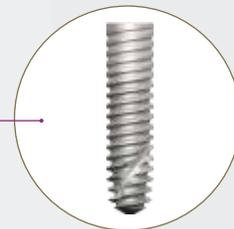

#### Hybrid contour

- Coronal: Cylindrical;
- Apex: Conical.

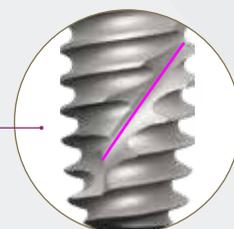


#### Active Apex

- Short tip;
- Helicoidal flutes.



### A SOLUTION FOR LIMITED BONE AVAILABILITY IN ALL BONE TYPES

Indicated for all bone types, the **Neodent® Helix GM Narrow** is specifically engineered to address esthetic challenges in situations with limited bone, thanks to its small diameter implant of 2.9mm.

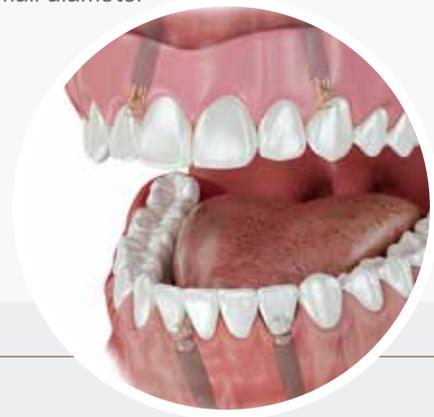

### COMPREHENSIVE PROSTHETIC PORTFOLIO FOR OPTIMIZED ESTHETIC AND FUNCTIONAL RESULTS


The **Helix GM Narrow** system was designed to offer clinicians greater levels of treatment flexibility with a comprehensive prosthetic portfolio, designed to meet patient expectations regarding short treatment times, esthetic and functional results.

It allows single and multi-unit restorations from screw and cement-retained, to removable prosthesis. The system also allows support for conventional and digital workflows supporting provide natural-looking restorations using either conventional or immediate protocols.



 <p>Titanium Temporary Abutment</p>	 <p>Titanium Base</p>	 <p>Universal Abutment</p>	 <p>Micro Abutment</p>	 <p>Attachment Removable</p>	 <p>Titanium Block (AG or Medentika Holder)</p>	
 Single-unit screw-retained prosthesis	 Single-unit cement-retained prosthesis	 Multiple-unit screw-retained prosthesis	 Multiple-unit cement-retained prosthesis	 Temporary	 Overdenture	 Straumann UNIQ workflow availability*

\*Check availability in your region

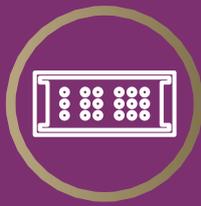


## Simplicity for treatment protocols

When it comes to implant procedures, both guided and non-guided surgeries can become unnecessarily complex due to the number of instruments that are usually needed. Additionally, the lack of clear instrument identification can increase the risk of errors and reduce patient satisfaction.

The **Helix GM Narrow** system provides an intuitive hybrid surgical kit designed to best suit any chosen surgical procedure, whether conventional or guided, adding even more simplicity to the system by using the Neo Screw connection.

### AN INTUITIVE AND FUNCTIONAL COMPACT SURGICAL CASSETTE

The **Helix GM Narrow** system allows intuitive conventional and guided surgeries with the functional compact surgical kit, to support improve outcomes and patient satisfaction.

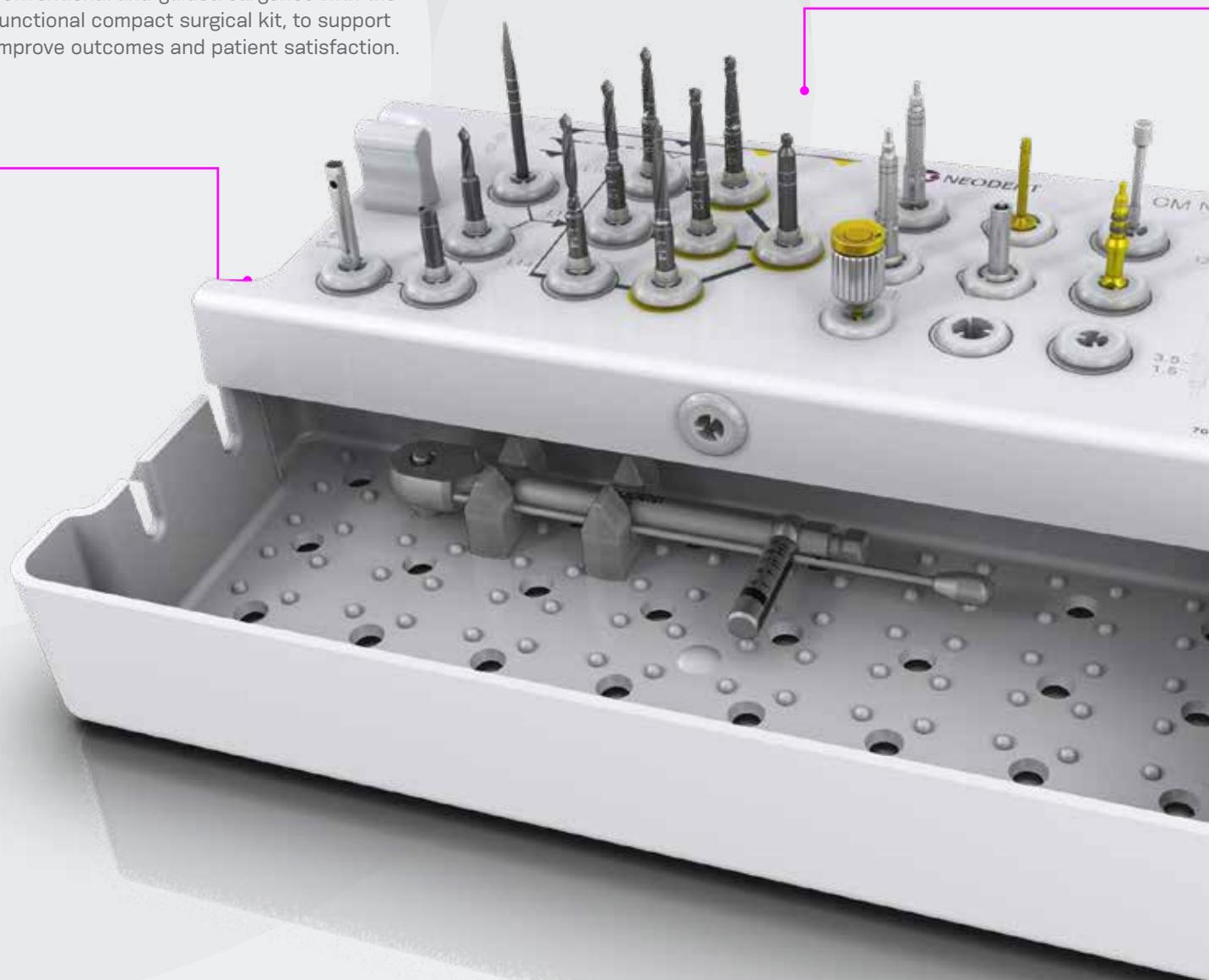



#### A PREDICTABLE GUIDED PROCEDURE WITH THE EASYGUIDE CONCEPT

The Neodent® EasyGuide concept offers straightforward guided surgery technique enabling surgical convenience with one-hand procedures, and pursuing predictable surgical results with confidence for accurate implant positioning.



#### ONE SCREWDRIVER AVAILABLE BOTH FOR NEODENT® GM AND GM NARROW

The **Helix GM Narrow** system features the Neo Screwdriver, which has a star attachment offering reliability and durability, compatible with all GM Narrow healing abutments and restorative screws.

# Clinical Case

**DR. GENINHO THOMÉ**


**BRAZIL**

Masters and PhD in Implantology  
Scientific Chairman Neodent®  
Chairman of the Board of  
Directors of Neodent®

**1** Initial clinical aspect



**2** Surgical aspect after implant and abutment placement

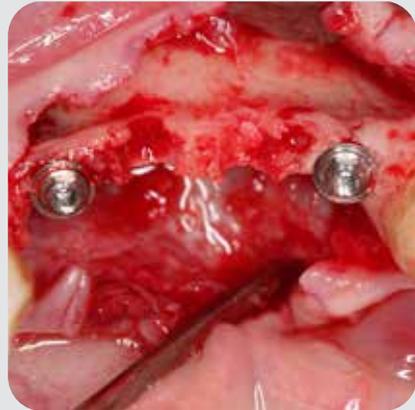

**3** Clinical aspect of surgical follow-up



**4** Radiographic aspect of surgical follow-up



**5** Final prosthetic rehabilitation



**6** Radiographic aspect of final prosthetic rehabilitation



# Smile through Life.

 **NEODENT.**  
A Straumann Group Brand

Neodent® is a global brand founded by a dentist for dentists, with the purpose of changing lives. Available in 95 countries, with a legacy of more than 30 years focused on ease of use, Neodent® Dental Implant Systems focus on progressive treatment concepts, such as immediacy with modern and reliable solutions to enable therapy access and affordability for creating new smiles every day.



## GLOBAL BRAND

Available in 97 countries, expanding our philosophy worldwide.



## FOUNDED BY A DENTIST FOR DENTISTS

A legacy of more than 30 years focused on ease of use.



## PROGRESSIVE TREATMENT CONCEPTS

Modern and reliable solutions.



## THERAPY ACCESS AND AFFORDABILITY

Accessibility to proven and affordable solutions.

- (1) SciELO - Brazil - Internal and Marginal Fit and Fracture Strength of Narrow Diameter Dental Implant-Abutment Assembly Internal and Marginal Fit and Fracture Strength of Narrow Diameter Dental Implant-Abutment Assembly
- (2) Al-Nsour MM, Chan HL, Wang HL. Effect of the platform-switching technique on preservation of peri-implant marginal bone: a systematic review. Int J Oral Maxillofac Implants. 2012 Jan-Feb;27(1):138-45.
- (3) Annibaldi S, Bignozzi I, Cristalli MP, et al. Peri-implant marginal bone level: a systematic review and meta-analysis of studies comparing platform switching versus conventionally restored implants. J Clin Periodontol. 2012 Nov;39(11):1097-113.
- (4) Hsu YT, Lin GH, Wang HL. Effects of Platform-Switching on Peri-implant Soft and Hard Tissue Outcomes: A Systematic Review and Meta-analysis. Int J Oral Maxillofac Implants. 2017;32(1):e9-e24.
- (5) Lazzara RJ, Porter SS. Platform switching: a new concept in implant dentistry for controlling postrestorative crestal bone levels. Int J Periodontics Restorative Dentistry. 2006 Feb;26(1):9-17.
- (6) Rocha S, Wagner W, Wiltfang J, Nicolau P, Moergel M, Messias A, Behrens E, Guerra F. Effect of platform switching on crestal bone levels around implants in the posterior mandible: 3 years results from a multicentre randomized clinical trial. J Clin Periodontol. 2016 Apr;43(4):374-82.
- (7) Novellino MM, Sesma N, Zanardi PR, Laganá DC. Resonance frequency analysis of dental implants placed at the posterior maxilla varying the surface treatment only: A randomized clinical trial. Clin Implant Dent Relat Res. 2017 Jun 20. doi: 10.1111/cid.12510. [Epub ahead of print]
- (8) Sartoretto SC, Alves AT, Resende RF, et al. Early osseointegration driven by the surface chemistry and wettability of dental implants. J Appl Oral Sci. 2015 May-Jun;23(3):279-87.
- (9) Sartoretto SC, Alves AT, Zarranz L, et al. Hydrophilic surface of Ti6Al4V-ELI alloy improves the early bone apposition of sheep tibia. Clin Oral Implants Res. 2016 Jun 17. doi: 10.1111/clr.12894. [Epub ahead of print]
- (10) Val JE, Gómez-Moreno G, Ruiz-Linares M, et al. Effects of Surface Treatment Modification and Implant Design in Implants Placed Crestal and Subcrestally Applying Delayed Loading Protocol. J Craniofac Surg. 2017 Mar;28(2):552-558.

JEDEN TAG EIN NEUES LÄCHELN

NUEVAS SONRISAS TODOS LOS DÍAS

NYA LEENDEN VARJE DAG

NEW SMILES EVERY DAY

NOVOS SORRISOS TODOS OS DIAS

CHAQUE JOUR DE NOUVEAUX SOURIRES

© 2025 – JIGC Indústria e Comércio de Materiais Dentários S.A. All rights reserved. All brands listed in this material, registered or not, are the property of JIGC Indústria e Comércio de Materiais Dentários S.A. The reproduction of this material without prior authorization is prohibited.

Neodent, Acqua, Helix GM, are trademarks or registered trademarks of JIGC Indústria e Comércio de Materiais Dentários S.A.

10797\_neodent\_ngm\_brochure\_en\_eU\_E00\_of\_300625

 **NEODENT**  
A Straumann Group Brand